


Sugita A V M Micro Clip

MIZUHO®

Catalog No.	Blade Length
-------------	--------------


07-932-02	2 mm
07-932-03	3 mm
07-932-04	4 mm
07-932-05	5 mm

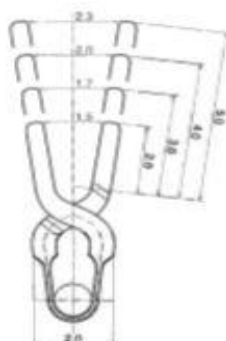

Clip can be easily picked up by an applicator while holding the container in one hand.



Sugita A V M Micro Clip Applicator

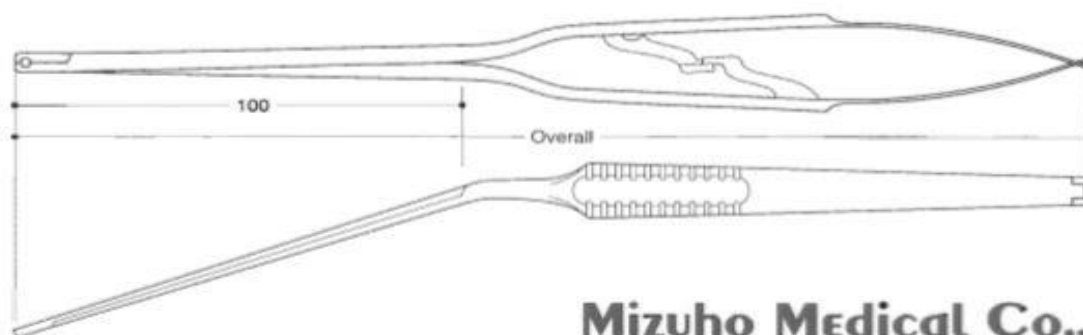

Catalog No.	Overall Length
-------------	----------------

07-935-10	21.2 cm
-----------	---------

Mizuho Medical Co., Ltd.

ICN building 5F 27-17,2-chome, Hongo, Bunkyo-ku,
Tokyo, 113-0033 Japan
Fax No.+81-3-3818-1705 Phone No.+81-3-3815-7101
E-Mail : marketing@mizuhomedical.co.jp

Mizuho America Inc.

133 Brimbal Avenue Beverly, MA 01915
Phone No.800-699-CLIP(2547) or 978-921-1718
Fax No.978-921-4003

Manufacturer

MIZUHO IKAKOGYO CO., LTD.
Tokyo, Japan



Sugita AVM Microclips are made of cobalt-chromium-molybdenum-based alloy having superior properties, specified under ASTM F1058 G2(ISO5832-7).

Being made of a non-ferromagnetic alloy, they are safe for MRI currently in clinical use.

To prevent damage and to facilitate handling, the clip is individually seated in a heat-resistant silicone sponge and packaged in a metal container, which is fully autoclavable.

Clips are made to have a holding force within the predetermined range of 50-70g.



Sugita AVM Micro Clip and Applier



Clip Packaging

