

CRANIOTOME

Pneumatic Bone Drilling / Cutting Device for Cranial Surgery

Features :

- Feel of Stable Gripping....
Grooves provided on the gripping section prevent inadvertent hand slip.
- Excellent Cutting Performance....
Two RPMs independently assigned for drilling and cutting devices afford unsurpassed high cutting performance.
Drilling Device 800 RPM · Cutting Device 25,000 RPM
- Double-Locking Mechanism....
Unlike conventional drilling / cutting devices of "insert-attachment" type, the double-locking mechanism has eliminated play of the bits or burs with less chances of breaking.
- A wide variety of cutting blades to choose from....

Specification :

- Pneumatic Supply Pressure · 8 kg/cm² (110 psi)
- Speed of Revolution · Drilling 800 RPM · Cutting 25,000 RPM
- Operating Gas · Nitrogen Gas (Air)





▲07-734-09
Pneumatic Bone Drill (Perforator)
Handpiece



▲07-734-01
Pneumatic Bone Cutter Handpiece



▲07-734-02
Bur Cover
(Fixed Bur Guard Type)



▲07-734-03
Bur Cover
(Without Bur Guard)



◀07-734-04
Pineapple Bur (Cylindrical Type)

◀07-734-05
Pineapple Bur (Conical Type)

◀07-734-07
Drill Bit for Preparatory Drilling (φ1.5)

◀07-734-06
Drill Bit for Preparatory Drilling (φ2.0)

◀07-734-08
Side Cutting Bur




▼07-734-10
Drill Bit for
Perforator
(Separable Type)



▼07-734-11
Inner Drill Bit (Disassembling Type)



▲01-460-00
Pressure Supply Hose

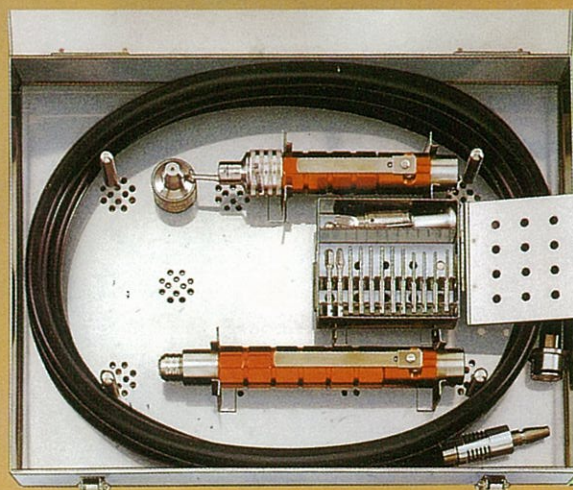

07-734-00 Contents of Set and Number of Pieces Included.

07-734-01	Pneumatic Bone Cutter Handpiece	1
07-734-02	Bur Cover (Fixed Bur Guard Type)	1
07-734-03	Bur Cover (Without Bur Guard)	1
07-734-04	Pineapple Bur (Cylindrical Type)	1
07-734-05	Pineapple Bur (Conical Type)	1
07-734-06	Drill Bit for Preparatory Drilling (φ2.0)	2
07-734-07	Drill Bit for Preparatory Drilling (φ1.5)	2
07-734-08	Side Cutting Bur	6
07-734-09	Pneumatic Bone Drill (Perforator) Handpiece	1
07-734-10	Drill Bit for Perforator (Separable Type)	1
07-734-11	Inner Drill Bit (Disassembling Type)	1
01-460-00	Pressure Supply Hose	1
07-734-50	Sterilizing Case	1

Option

01-461-00

Gas Regulator

07-734-20

Bur Guard, Rotatable Type Bur Cover



Mizuho Medical Co., Ltd.

ICN Building 5F 27-17, 2-chome, Hongo, Bunkyo-ku, Tokyo
113-0033 JAPAN

Phone No. 03-3815-7101 Fax. No. 03-3818-1705

Manufacturer

MIZUHO IKAKOGYO CO., LTD.

Tokyo, Japan

2003/02/2000/1